



PCR WORK STATION



Manufacturer of
Life Science Research
Equipments



EPS offers PCR workstation which has self-contained work area that will help protect your PCR runs against contamination. UV bulbs irradiate the work area prior to use, reducing the possibility that contaminating DNA will be amplified. Containment features reduce the chance of airborne contamination. The non-ventilated hood has acrylic windows for a clear view, an acid-resistant SS work surface.

PCR amplification is extremely sensitive to contamination, prevention of contamination requires good laboratory practices to minimize external or cross-contamination during reagent and sample preparation.

Applications

- Medicine
- Environmental Study
- Crime Scene Investigation

Construction

- MS Power coated with Stainless Steel Working Top
- Door: PC Sheet for UV Resistance
- Power Socket Optional



Technical Specifications & Ordering Information:

MODEL	EPS/PCR/600W	EPS/PCR/900W	EPS/PCR/1200W
Chamber size (feet) LxDxH	2'X2'X2.5'	3'X2'X2.5'	4'X2'X2.5'
Cabinet	1X15 watts	1X15 Watts	2X15 Watts
UV source	1X15 watts	1X15 Watts	2X15 Watts
Wavelength	254 nm	254 nm	254 nm
Florescent lamp	1X8 Watts	2X8 Watts	2X8 Watts
Pipet Holder	5 placings	5 placings	5 placings
Timer	Optional	Optional	Optional
Power	230V AC, 50Hz	230V AC, 50Hz	230V AC, 50Hz



UDYAM REGISTRATION CERTIFICATE	
UDYAM REGISTRATION NUMBER	UDYAM-TN-42-613876
NAME OF ENTERPRISE	MS EPS BIOSOLUTIONS
TYPE OF ENTERPRISE *	MICRO (Based on FY 2019-20)
MAJOR ACTIVITY	MANUFACTURING
SOCIAL CATEGORY OF ENTREPRENEUR	OBC
NAME OF UNIT(S)	1. EPS BIOSOLUTIONS
OFFICIAL ADDRESS OF ENTERPRISE	Flat No. 4, 1st Floor, 1st Cross Street, Chittibabu Nagar, 1st Cross Street, Chittibabu Nagar, Chennai, Chengalpattu, Tamil Nadu, 600100
DATE OF INCORPORATION / REGISTRATION OF ENTERPRISE	16/11/2008
DATE OF COMMENCEMENT OF PRODUCTION/BUSINESS	16/11/2008
NATIONAL INDUSTRY CLASSIFICATION CODE(S)	31.20 - Other manufacturing
DATE OF UDYAM REGISTRATION	06/06/2022

* In case of graduation (upward movement) of status of an enterprise, the benefit of the Government Scheme will be available on the basis of the date of graduation.

For any assistance, you may contact:

1. District Industries Centres: CHENNAI (TAMIL NADU)

2. MSME DE: CHENNAI (TAMIL NADU)


BE A CHAMPION with the Ministry of MSME

Certificate of Registration of Trade Mark, Section 23 (2), Rule 56 (1)	
Trade Mark No.	3938903
Class	9
Under No.	3938903
in of the date	06 September 2018
in respect of	Goods and service as annexed
Trade Mark as annexed	
Registered on	04th day of April, 2019
Registrar of Trademarks	MUMBAI

OUR CLIENTS



Our Dealers Network

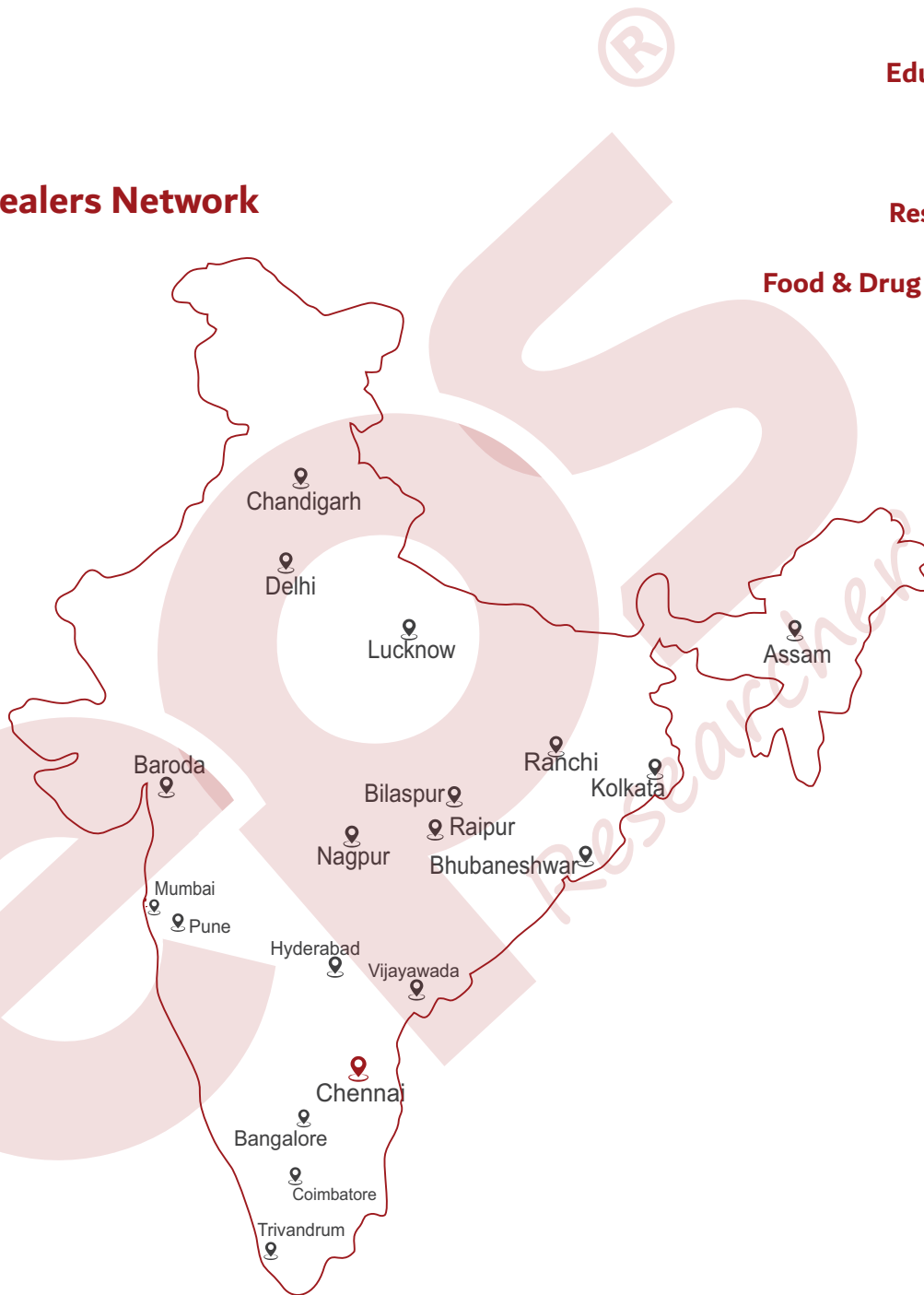

Educational Institutions ◀

Pharma Industries ◀

Research Organizations ◀

Food & Drug Testing Laboratories ◀

Diagnostic Labs ◀

Aquaculture Labs ◀

EPS Biosolutions

Registered Office:

Plot 4B, Annai Abode Apartment, 1st Cross Street,
Chitti Babu Nagar, Pallikaranai, Chennai - 600 100.

Tamilnadu. India. ☎: 044-4207 2939.

🌐 info@epsbiosolutions.com | sales@epsbiosolutions.com

✉ www.epsbiosolutions.com

Factory:

No. 28/17, Bajanai Kovil Street,
Velachery, Chennai - 600 042. Tamil Nadu, India

Off - Orandiamman Kovil Street,

☎ 044-42072939 🌐 info@epsbiosolutions.com

✉ www.epsbiosolutions.com